





ELECTROPHORUS – handle & plates

EM1780-001 (50mm diam) EM1781-001 (75mm diam) EM1782-001 (150mm diam)


DESCRIPTION:

The classic Electrophorus is a metal disc with insulating handle, a metal base plate and an insulating plate that rests on the base plate. When an electric charge is applied between the base plate and the disc, a capacitor is formed and the plate can be removed with an electric charge residing.

EM1780-001 50mmD EM1781-001 75mm D

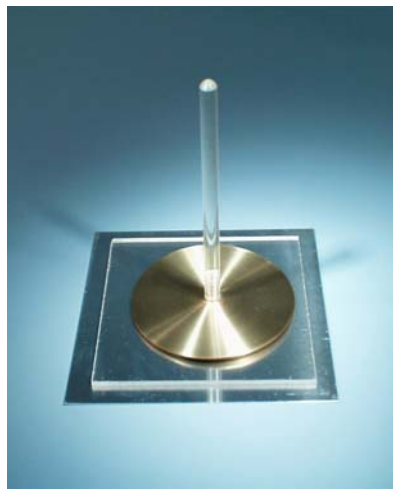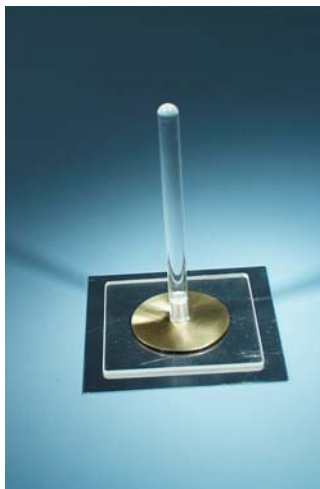

EM1782-001 150mm D

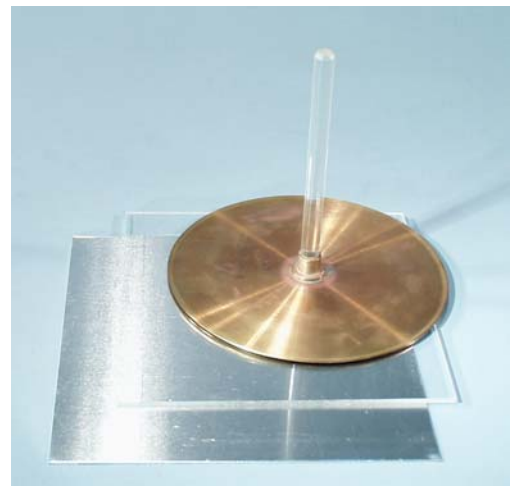

Physical sizes:

EM1780-001	Disc: 50mm diameter	Weight: 0.1 kg
EM1781-001	Disc: 75mm diameter	Weight: 0.15 kg
EM1782-001	Disc: 100mm diameter	Weight: 0.2 kg

Designed and manufactured in Australia