

POLAROID FILTERS – pair with scale

Cat: HL1883-050 (100mm sq. with circular scale)


DESCRIPTION:

The standard IEC Linear Polarising films are available in pairs of squares of various sizes. They can be used as squares or they can be cut into shapes or strips to suit various experiments.

The transmission is 38% when polaroids are aligned. When misaligned at 90°, the transmission reduces to less than 0.2% (almost complete blackout).

This pair of polarised filters is 100mm square and each is printed with a circular scale to make the angular alignment between them easier to understand.

HL1883-050 polaroid filters with scale

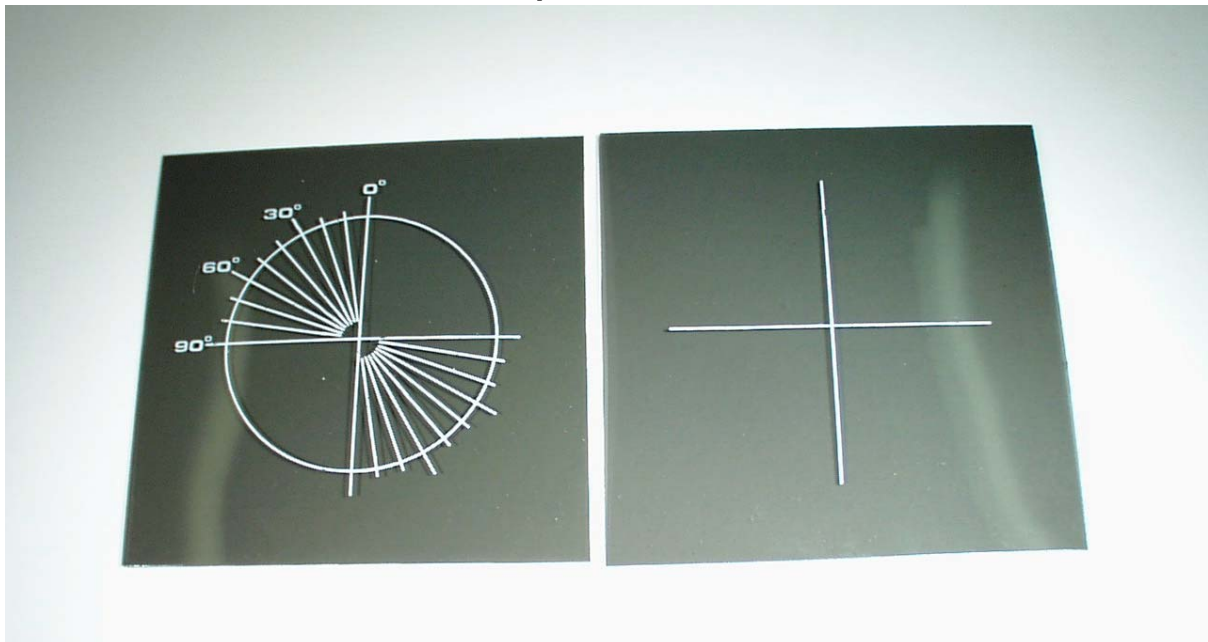

Code	Type	Size & thickness
HL1883-050	Linear polariser	Pair of squares with scale 100x100mm x 0.18mm thick

Designed and manufactured in Australia